

Article No. S40320



Preparation Products

# Strip Cutter CUT7

Special double knife cutter for thinnest papers and foils.

For:




✓ PAPER



✓ FOILS



✓ TISSUE

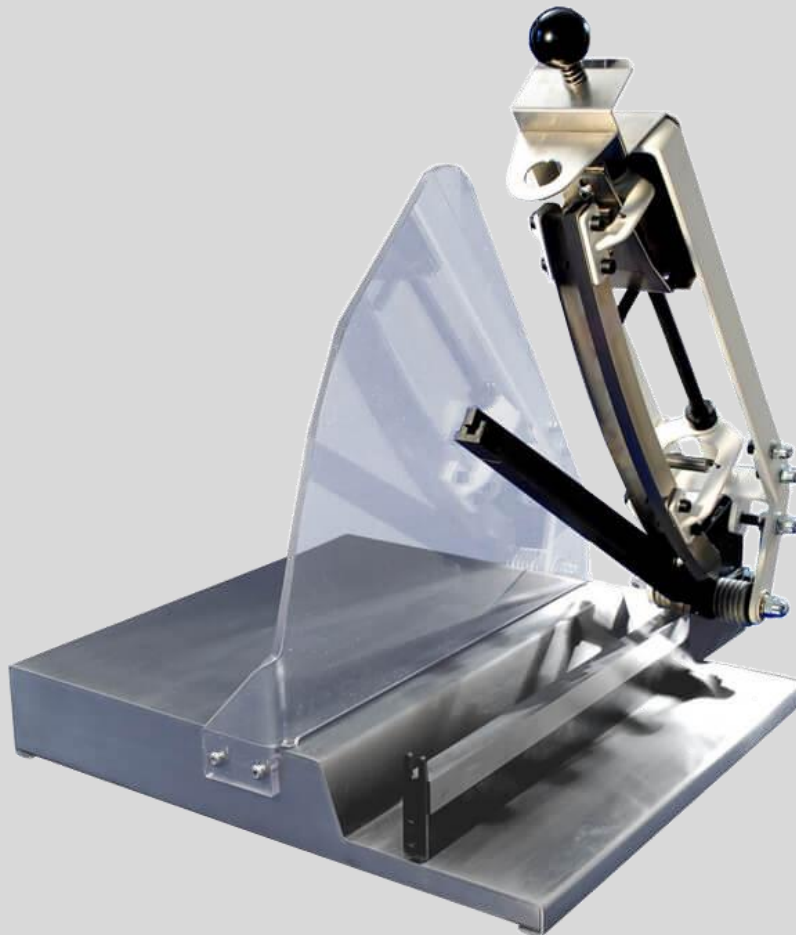

## MOST IMPORTANT BENEFITS:

- ✓ For thinnest paper and foils
- ✓ Cutting lock – prevents injuries and fixes the material
- ✓ Rounded blades ensure a precise cut
- ✓ Parallel guided blades ensure highest parallelism of the cut samples

 **FRANK-PTI**  
QUALITY TESTING INSTRUMENTS

🏠 FRANK-PTI GmbH  
Auf der Aue 1  
69488 Birkenau

☎ +49 (0) 6201 84 0  
📠 +49 (0) 6207 84 290  
✉ office@frank-pti.com

🌐 www.frank-pti.com

## PRODUCT DESCRIPTION

Special double knife cutter for thinnest papers and foils.

The Cut7 cuts paper from 8 g / m<sup>2</sup> and is used for the thinnest papers (8-30g / m<sup>2</sup>). Even the thinnest foils and aluminum foils can be cut with it. A cutting lock prevents injuries and additionally clamps the paper so that a straight cut is made and the material does not slip.

When ordering please specify the material to be cut, as the cutter is specially adjusted for this.



Spring-loaded round blades for a precise and clean cut

## TECHNICAL DATA

### DEVICE / INSTRUMENT

- Suitable for the finest papers and foils
- Cutting lock prevents injuries and fixes the material to be cut
- Spring-loaded rounded blades for a precise, clean cut
- Parallel guided blades ensure highest parallelism of the cut samples
- For cut material between 8-30g / m<sup>2</sup>
- Specially developed for flat materials:  
Aluminum, copper, brass foil  
Paper, nonwovens, textile and similar materials
- Cam plate guided hold-down device
- Cutting length max. 300mm



### DIMENSIONS

S403200000	<b>L x H x B</b> 350 x 350 x 550 mm net
Weight:	approx. 12 kg

### ARTICLES / MODELS

S403200000	Cutting length 15,0 mm
S40320000x	Cutting length 50,0 mm
	Cutting length max. 300mm

\*more models, individual sizes and materials on request

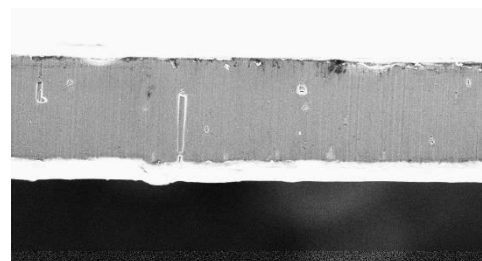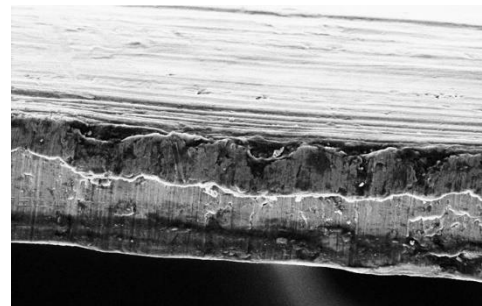

Comparison: regular strip cutter above / CUT7 strip cutter below